

# Cloud-based Signature Application

Mobile Access to the SOMS and Signature Planning Software

Powered by:



## Overview

With the use of this cloud-based tool, users can now access the Signature Online Management System (SOMS) and Signature Planning Software, using virtually any platform or device. Platforms and devices include:

- Macintosh® (OS)
- PC (Windows®)
- iPad®
- Android® Devices



## Getting Started

Complete the attached registration form attached here (or on mybiomet) and submit to **signature@biomet.com**. You will receive a Username and Password. All EMEA users must use the following **<Host URL>**:

Europe: **<Host URL>** = BiometEU.MyVpresence.com

## Platforms and Devices

Please refer to the appropriate section below for your specific device or operating system.


### PC (Microsoft® Windows®)

1. Using your web browser, go to **http://<Host URL>**.
2. The first time you visit this page, you will be prompted to install a plug-in. Be sure to select the **"Allow Plug-ins"** prompt if it is displayed. There will be a longer than typical delay while the plug-in is downloaded. After the plug-in is loaded, you will be prompted for your **Username** and **Password**.
3. After entering your username and password, the GO-Global application window will appear with the Internet Explorer icon (SOMS) and the Signature System icon (Planning Software).
4. Double-click the appropriate icon to launch the desired application.

### Macintosh® (Apple® OS X)

1. Using your web browser, go to **http://<Host URL>/clients.html**.
2. Select the **Mac OS X Client** option that is presented.
3. Follow the on-screen instructions to download and install the plug-in.
4. After the plug-in is installed, open your web browser and go to **http://<Host URL>**, where you will be prompted for your **Username** and **Password**.

**Note:** After installation, the GO-Global icon can be dragged from the Applications folder into the desktop dock. To launch the application, click the GO-Global icon, input the **<Host URL>** and click "Connect."

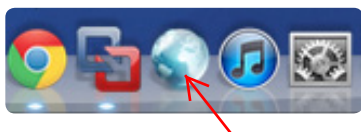


## Macintosh® (Apple® OS X) cont.

5. After entering your username and password, the GO-Global application window will appear with the Internet Explorer icon (SOMS) and the Signature System icon (Planning Software).
6. Double-click the appropriate icon to launch the desired application.


**Note:** If using the trackpad on a MacBook, a two-finger-click moved left and right on the trackpad, will rotate the view where the cursor is located. Sliding two fingers up and down on the trackpad (without clicking) will zoom in or out on the cursor location.

## Apple® iPad®

1. This cloud-based tool uses the GO-Global app (  ) on the Apple iPad. Download the GO-Global app from the Apple App Store at <http://itunes.apple.com/us/app/go-global/id441263366?mt=8>.
2. When running the GO-Global app for the first time, click the plus icon next to the “Host List.”




The “Add Host” pop box will appear. Select the “Name” button next to the “Host Address Type” and type in the **<Host URL>** web address next to the “Host Name.” Next, enter your **Username** and **Password** to configure a host connection.

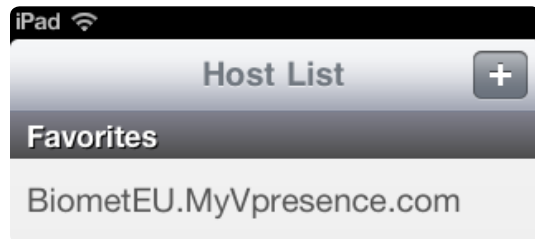

**Note:** Once configured, this **<Host URL>** will be listed in your Favorites. On each subsequent visit, select the saved URL listed in your Favorites followed by the “Connect” button in the upper right-hand corner.

3. After entering your username and password, the GO-Global application window will appear with the Internet Explorer icon (SOMS) and the Signature System icon (Planning Software).
4. Double-click the appropriate icon to launch the desired application.


You can also access the GO-Global User Guide by clicking the information icon (  ), which is located at the bottom of the Host List page.

### Common Finger Swipes:


- Enable Keyboard - Three finger simultaneous tap on display
- Rotate View - Use two fingers to click on the view, pause slightly, and then drag fingers left and right
- Zoom - Use two fingers to click on the view and immediately drag fingers up or down without pausing or delay
- Dragging Features – Use a single finger to select an item and drag it left or right

## Android® OS

1. This cloud-based tool uses the GO-Global app (  ) on Android® devices. Download the GO-Global app from the Android® app market at: **<https://market.android.com/details?id=com.graphon.goglobal>**.
2. When running the GO-Global app for the first time, click the plus icon next to the “Host List.”

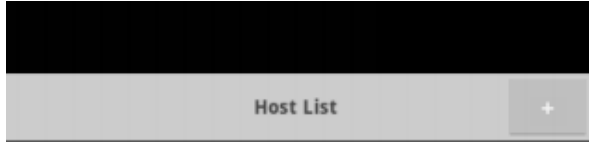


Use the **<Host URL>** web address and your **Username** and **Password** to configure a host connection.

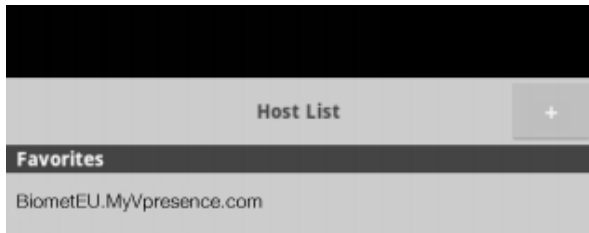

**Note:** Once configured, this **<Host URL>** will be listed in your Favorites. On each subsequent visit, select the saved URL listed in your Favorites followed by the “Connect” button in the upper right-hand corner.

3. After entering your username and password, the GO-Global application window will appear with the Internet Explorer icon (SOMS) and the Signature System icon (Planning Software).
4. Double-click the appropriate icon to launch the desired application.

You can also access the GO-Global User Guide by clicking the information icon (  ), which is located at the bottom of the Host List page.

### Common Finger Swipes:

- Enable Keyboard - Three finger simultaneous tap on display
- Rotate View - Use two fingers to click on the view, pause slightly, and then drag fingers left and right
- Zoom - Use two fingers to click on the view and immediately drag fingers up or down without pausing or delay
- Dragging Features – Use a single finger to select an item and drag it left or right

## Support

For additional questions or support, please contact the Signature Customer Support Team at:

- Phone: **+1 574.371.3710**
- Email: **[signature@biomet.com](mailto:signature@biomet.com)**

Apple and iPad are registered trademarks of Apple, Inc.  
Android is a registered trademark of Google, Inc.  
Virtuon is registered trademark of Virtuon, Inc.

CONFIDENTIAL material intended solely for Biomet employees and the sales force.

### **European Central Marketing**

Waterton Industrial Estate  
Bridgend  
South Wales  
CF31 3XA  
+44 [0] 1656 655221  
+44 [0] 1656 645454

### **Distributed in the UK by:**

Biomet UK  
Waterton Industrial Estate  
Bridgend  
South Wales  
CF31 3XA  
©2012 Biomet Orthopedics

**[www.biomet.com](http://www.biomet.com)**



**One Surgeon. One Patient.**

# Cloud-based Signature Application

Mobile Access to the SOMS and Signature Planning Software

Powered by:



## User Registration Form

To receive your Username and Password for the Cloud-based Signature Application, please email this completed form to **signature@biomet.com**. If you have any questions, please contact the Signature Customer Support Team at **+1 574.371.3710**.

Name:

User Type:

Distributor:

Country:

Email:

Phone:

