



X - PSI Image Acquisition *


Position X-Ray Calibration Strap and Marker 3D

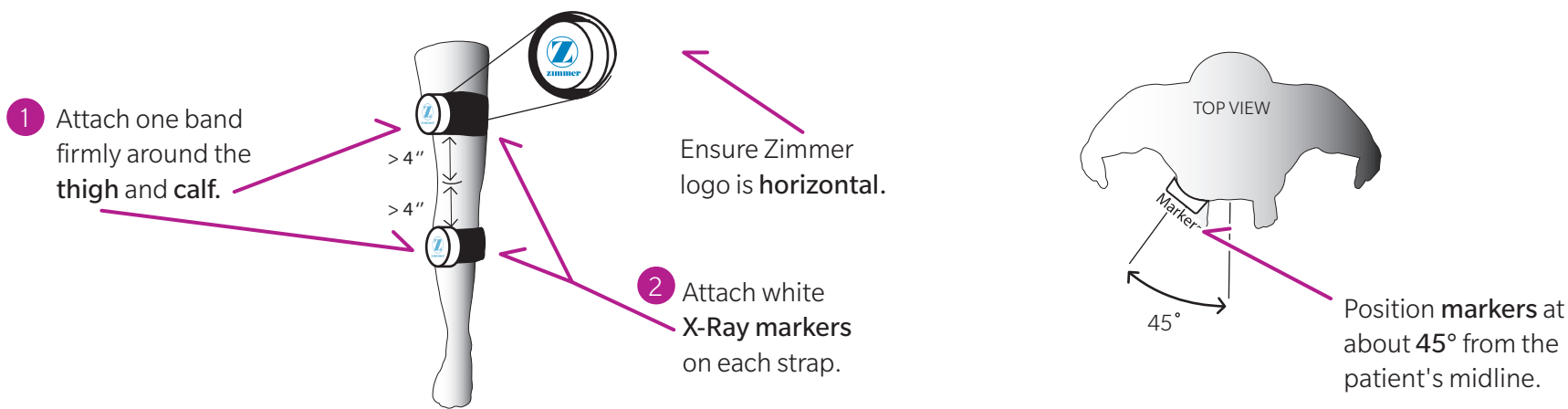

Position patient and take A/P full leg

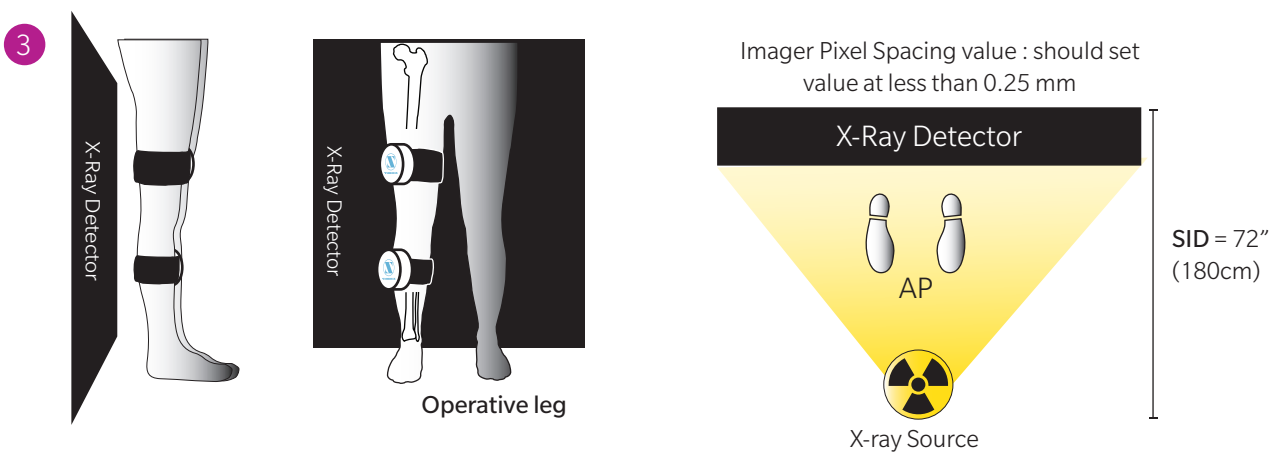

Position patient and take LAT full leg

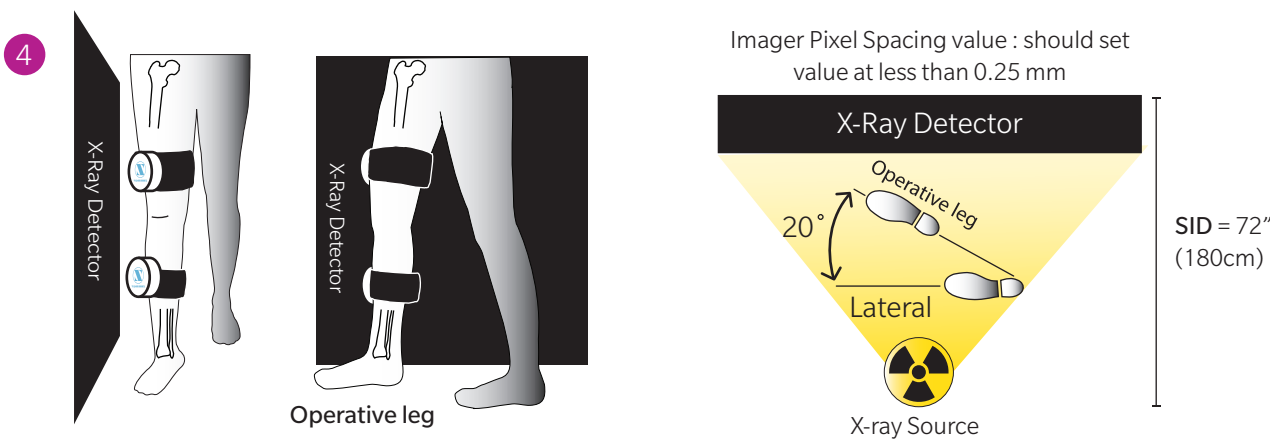

5 TIPS:

- In standing position, capture the entire leg image from above the femoral head to below the ankle joint .
- Prevent patient movement between any of sequential images in AP or LAT.
- No repositioning of the Markers is permitted during the procedure.
- SID (Source to Image Distance) value: must be fixed during the entire study.


Prepare and transfer final images to Zimmer Biomet

- 6 The following should be available on the DICOM tag or engraved on images:

SID
Imager Pixel Spacing
Patient Gender
Laterality
PSI Case ID
Surgeon Name

- 7 Transfer DICOM Images (final stitched and non-stitched images) to Zimmer

Sample of good final stitched images

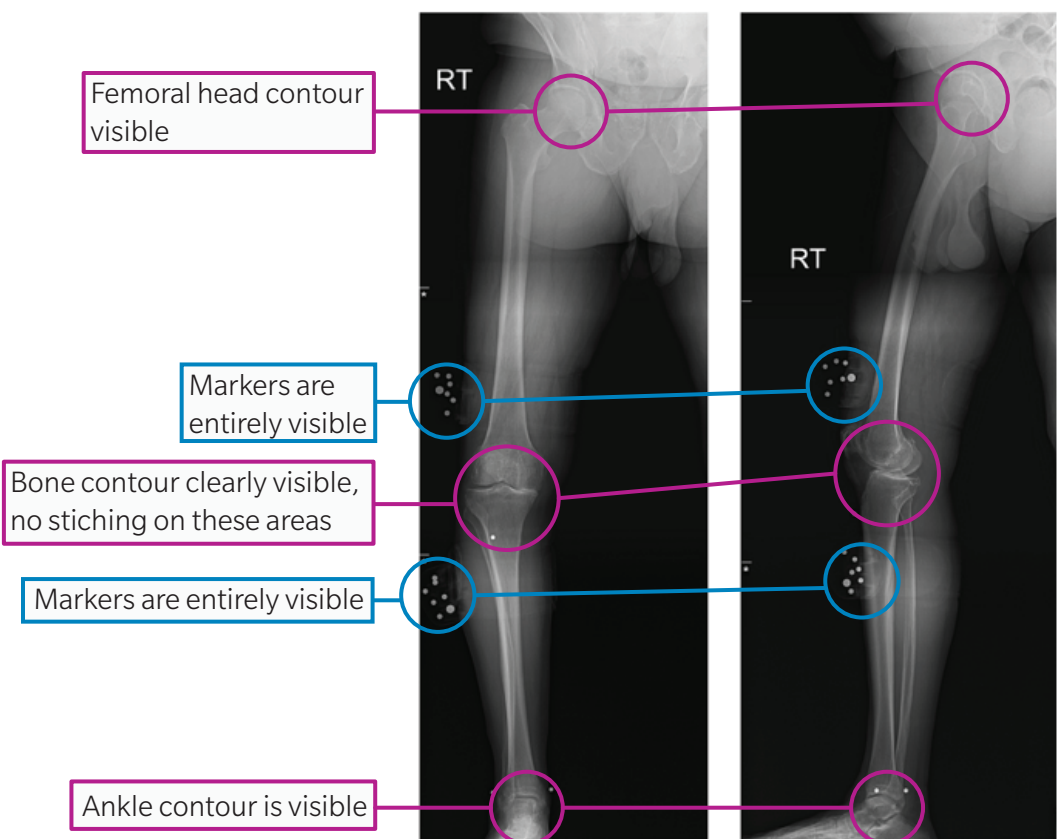

*For more details on each step please refer to XPSI Acquisition protocol (803.305)